


A-dec 300 Light^{LED}



 a dec[®]

Pure brilliance.

The first thing you'll notice about A-dec LED dental lights is the pure, neutral white light. Add our touchpad and all of the light modes and functions are available from one place. Your entire operator is bright, clutter free and connected, which makes it easy on your eyes—and your patient's.

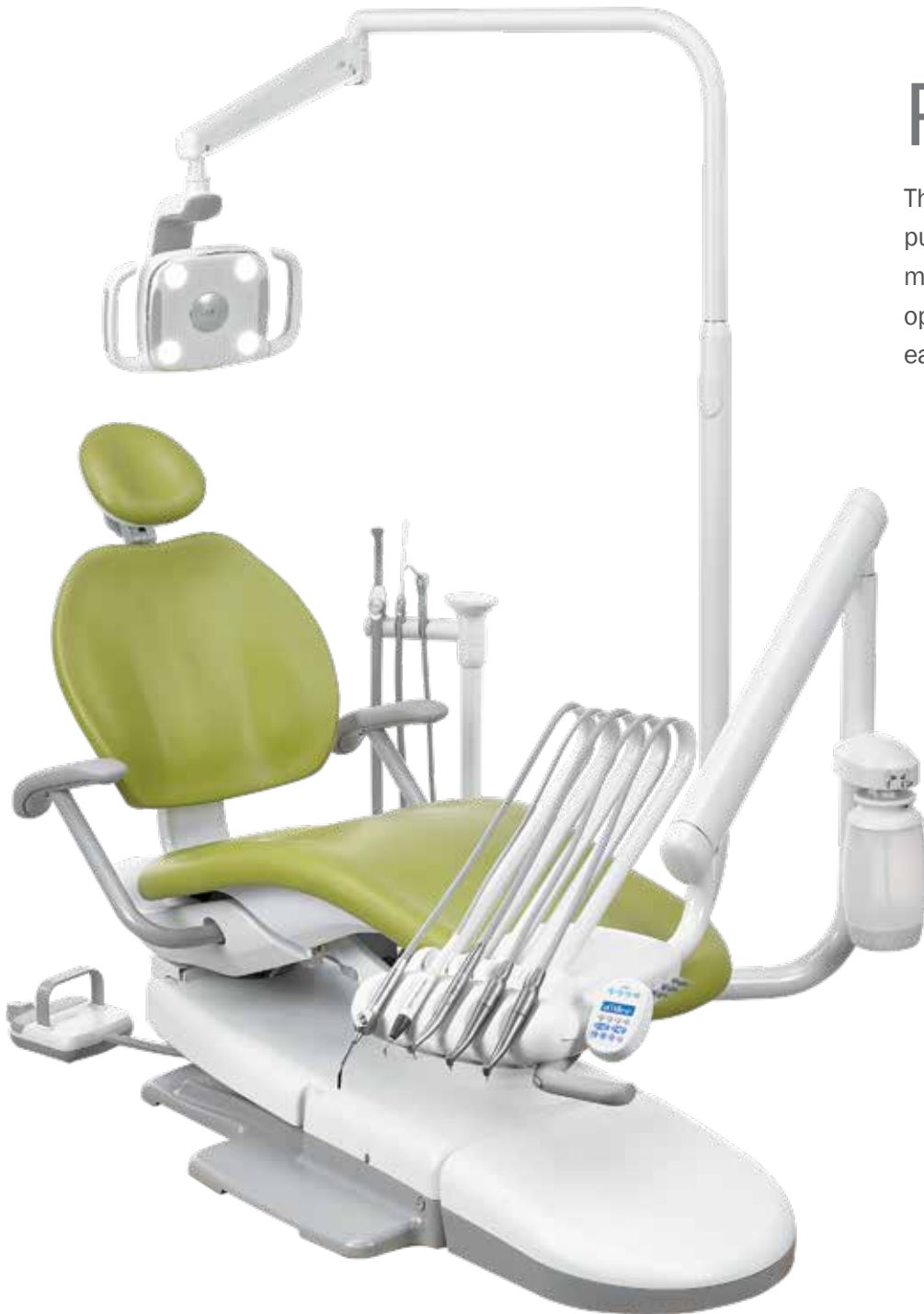

Backed by a 5-year warranty.

The difference is in the details.



PURE INNOVATION

Designed with silent passive cooling so the light stays cool—and so do you.



SMART DESIGN

Designing and manufacturing in house means we control quality. That's why the large, ergonomic handles aren't just easy to grip: they're made from industrial grade resin that resists impact damage.



UNPARALLELED TESTING

"Good enough" just isn't. A-dec designs and tests our lights for a minimum 20-year life. So the light always works when you need it.

A-dec LED


What color is natural?

Seeing colors as nature intended means an accurate representation and mix of the entire color spectrum. With true-to-life tones, a clear diagnosis is that much easier.

Solid performance. A-dec lights begin with proven parts, optically engineered specifically for dentists in their environment. A-dec LED lights are backed by a five-year warranty, but designed for a minimum 20-year life.

Easy on the eyes. Innovative light distribution creates a uniform pattern with a “stadium” effect to minimize shadowing, and feathered edge light pattern that reduces eye fatigue.

Exceptional illumination. A high color rendering index (CRI) of 94 mimics sunlight’s clarity and floods the oral cavity with light that reflects colors accurately for soft and hard tissue diagnosis. This high-quality light pattern conforms to the latest ISO9680 standards for brightness, color, rendering and uniformity.

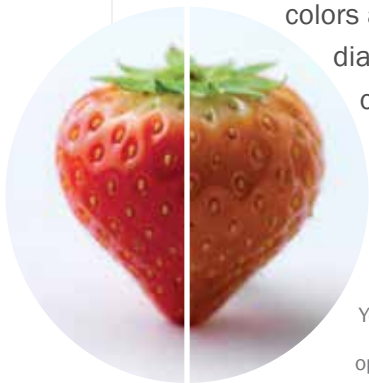

High CRI

Low CRI

You can see the difference when compared to less refined LED operator lights with a lower color rendering index (CRI).





Low power. No noise. No noise is good noise. The A-dec 300 uses a unique passive cooling method so there's no fan, and no radiant heat.

Ecologically and economically smart. A-dec LEDs consume 85% less power than halogen, without sacrificing brightness. Low heat and free of toxic chemicals, it produces little infrared light and close to no UV emissions. And, with a life expectancy of 40,000 hours, you'll never have to change a bulb.

Perfect white light. Neutral is everything. The color of light, measured in Kelvins, affects everything you do as a doctor. Consistent 5,000K illumination produces a balanced clean, white light that replicates natural sunlight.





NEW A-dec 300 LED

There is a light that's just right for everyone. The new A-dec 300 was designed specifically for those who want the clarity and accuracy of the A-dec 500 LED, in a simple, economical version that's anything but a compromise.

A-dec 300 LED Light

- A-dec 300 LED light choices
- Integrated (data communication system “DCS”) with A-dec touchpad. *Controlled via push button on light head or through integration with A-dec touchpad.*
 - Non-integrated with manual switch. *Limited to control via push button on light head. Available only on Performer or A-dec 200.*

| Feature | Description | Why it's important |
|-----------------------------|---|--|
| Stadium lighting | 4 bright white LEDs | Reduces shadowing |
| Light intensity (two modes) | High: 25,000 lux (2323 fc) Low: 8,000 lux (743 fc) | Enough light to do the job right |
| Thermal output | 51 BTUs (halogen average 325 BTUs) | No radiant heat |
| Composite mode | Yes - 8,000 lux (slow cure) | Allows more work time |
| Color temperature | 5,000K | Accurate color matching |
| Light pattern | 3.8" x 5.7" at 27.6" focal distance (95 mm x 145 mm at 700 mm) | Lights oral cavity, not patient's eyes |
| Feathered edge | Yes | Reduces eye fatigue |
| Lifespan | 40,000 hours (~20 years) | Never change a bulb |
| Rotational axes | 3-axis: horizontal, vertical, diagonal | Choose optimal ergonomic positioning |
| Compatibility | A-dec 500, A-dec 400, A-dec 300, A-dec 200, Performer | Select the option that works with your equipment |

Note: Composite mode is an industry accepted term, not intended to meet the 2014 ISO 9680 5.2.10 standard.



A-dec LEDs provide uninterrupted performance for 20 years, while consuming 75-85% less energy than halogen technology.

See the advantages.



UNLIMITED POSITIONING

540° of horizontal rotation mean that you can rotate the light 1½ times in either direction. Combine that with 120° head tilt and 80° diagonal adjustment, and perfect placement is a sure thing.



TOUCHPAD CONTROLLED

Programmable buttons automatically turn light on when the chair back reclines, and powers it off when the chair moves to entry/exit position (*requires A-dec touchpad*).



EASY ON/OFF

A simple bump lets you control the light, even during procedures.



CLEAN AND SIMPLE

With asepsis in mind, the light shield fits flat and tight to prevent dust accumulation inside. Barrier ready touch surfaces make it easy to clean the outside, too.



Light Options

What works for you? Choose a mounting option that maximizes your working style.

All options may not be available.
Check with your dealer for availability.



Track mount



Ceiling mount



Radius® mount



Chair mount



Wall/cabinet mount



A-dec Headquarters

2601 Crestview Drive
Newberg, Oregon 97132 USA
Tel: 1.800.547.1883 within USA/CAN
Tel: +1.503.538.7478 outside USA/CAN
Fax: +1.503.538.0276
a-dec.com

A-dec Australia

Tel: 1.800.225.010 within AUS
Tel: +61.(0)2.8332.4000 outside AUS

A-dec China

Tel: 400.600.5434 within CN
Tel: +86.571.89026088 outside CN

A-dec United Kingdom

Tel: 0800.233285 within UK
Tel: +44.(0).24.7635.0901 outside UK

A-dec, the A-dec logo, A-dec Inspire, Cascade, Century Plus, Continental, Decade, ICX, ICV, Performer, Preference, Preference Collection, Preference ICC, and Radius are trademarks of A-dec Inc. and are registered in the U.S. and other countries. A-dec 500, A-dec 400, A-dec 300, A-dec 200, and reliablecreativesolutions are also trademarks of A-dec Inc. None of the trademarks or trade names in this document may be reproduced, copied, or manipulated in any manner without the express, written approval of the trademark owner.

Certain touchpad symbols are proprietary to A-dec Inc. Any use of these symbols, in whole or in part, without the express written consent of A-dec Inc. is strictly prohibited.

The colors and products shown may vary from the actual materials and are subject to change without notice. Please see your authorized A-dec dealer for accurate color samples and the most current product information.